

## Deltalab Natural Polyethylene Caps for 15 - 16mm Tubes (with Internal Chamber)

**Product Code:** 305002


**Brand:** Deltalab

Ideal for capping and recapping 15 - 16mm blood collection tubes, Deltalab's Natural Polyethylene Caps for 15 - 16mm Tubes (with Internal Chamber) are designed to be placed and removed securely. Each cap has an internal chamber for more liquid filling, thus compensating for the loss of space inside the tube caused by the insertion of a cap. They are ergonomical and have a ribbed surface. Their shape is slightly conical for better handling.

These leakproof caps also feature an external flap to avoid contact between liquid and fingers when opening the tube, prevent dripping when the cap is left on the bench top, and provide easy single-handed opening.

### Specifications:

<b>Colour</b>	Natural
<b>Internal Diameter (mm)</b>	13.95
<b>External Diameter (mm)</b>	19.30
<b>Total Length (mm)</b>	13.15
<b>Material</b>	Polyethylene
<b>Sterility</b>	Non-Sterile
<b>Quantity</b>	2500

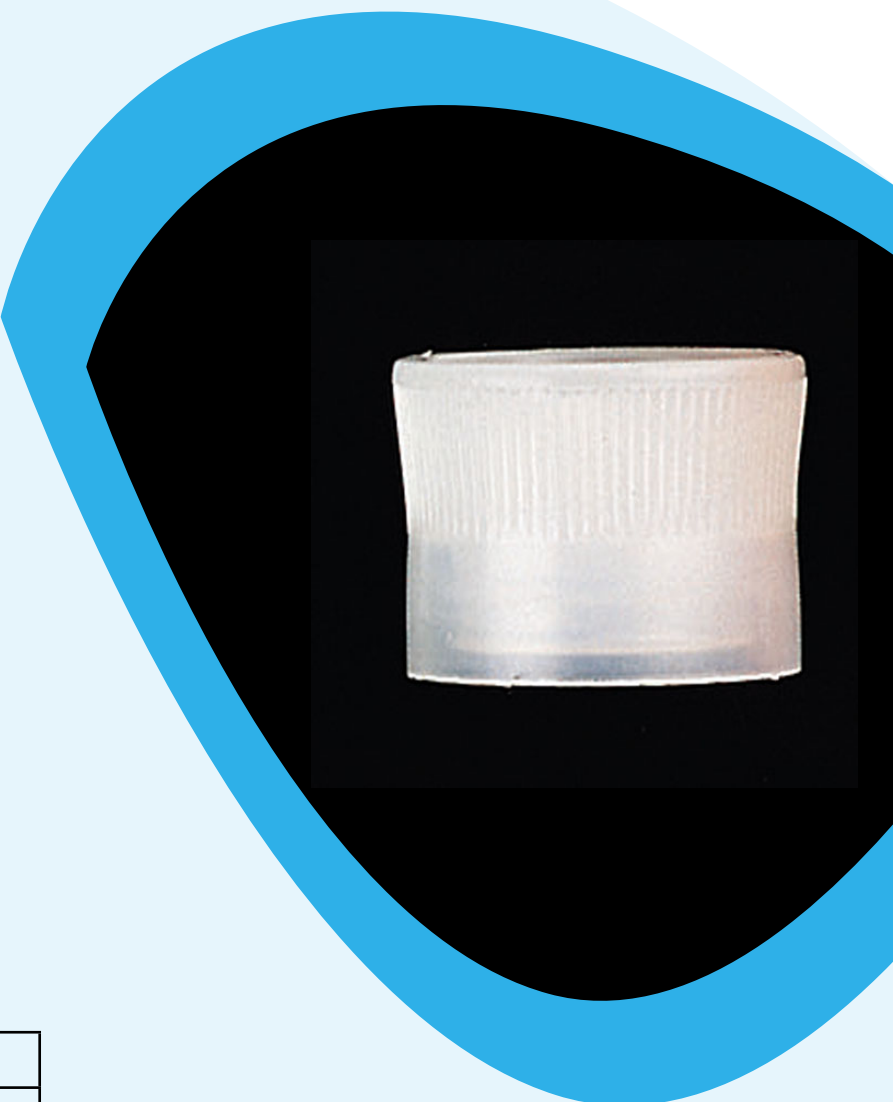

Unit 3, Tower Business Park  
Warpsgrove Lane  
Chalgrove  
Oxfordshire  
OX44 7XZ

**Tel:** (+44) 01865 400321

**Email:** [enquiries@medlinescientific.com](mailto:enquiries@medlinescientific.com)

**Website:** [www.medlinescientific.com](http://www.medlinescientific.com)