


## Deltalab 120ml Containers for Vacuum Systems (In Bulks)

**Product Code:** 408702.G

**Brand:** Deltalab

Designed to be used with vacuum tubes, Deltalab's 120ml Containers for Vacuum Systems (In Bulks) allow a clean transfer of the urine from the container to the vacuum tube. Manufactured in transparent polypropylene and with leakproof caps made of yellow polyethylene, these laboratory containers feature an internal obturating ring. They also include a cannula with a needle inside protected by a rubber cap and a plastic suction nozzle. The caps feature a label that seals the cannula in order to prevent contact with the collection needle. The containers also feature a wide marking by the cannula.

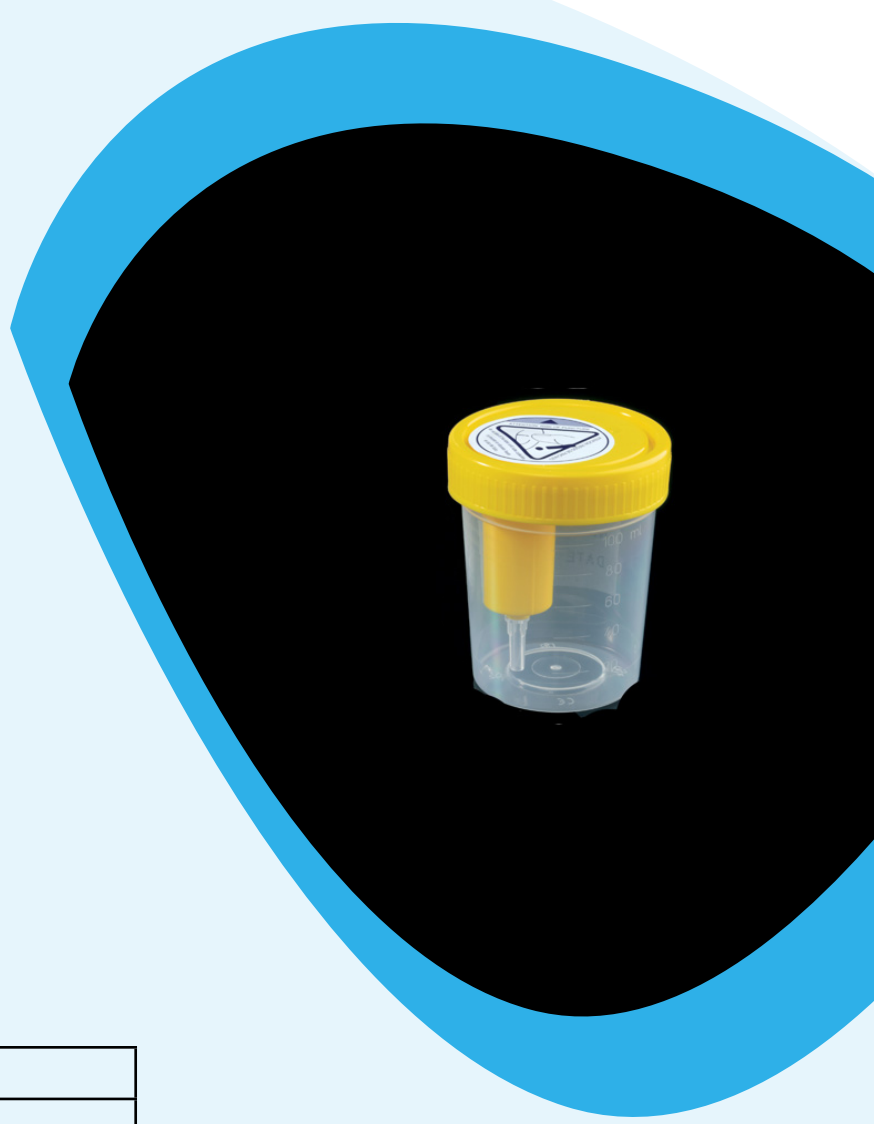

### Specifications

<b>Max. Volume (ml, with cap)</b>	120
<b>Recommended Volume (ml)</b>	60
<b>Dimensions (mm)</b>	57 x 73
<b>Presentation</b>	In bulk
<b>Sterility</b>	Sterile
<b>Quantity</b>	250

Unit 3, Tower Business Park  
Warpsgrove Lane  
Chalgrove  
Oxfordshire  
OX44 7XZ

**Tel:** (+44) 01865 400321

**Email:** [enquiries@medlinescientific.com](mailto:enquiries@medlinescientific.com)

**Website:** [www.medlinescientific.com](http://www.medlinescientific.com)