

Deltalab 2ml Non-Skirted Sterile Cryovials - With Internal Threads

Product Code: 409002

Brand: Deltalab

Made of autoclavable polypropylene, Deltalab's 2ml Non-Skirted Sterile Cryovials - With Internal Threads are designed for the storage of biological material at temperatures as low as -196°C (liquid nitrogen but in gas state).

Each cryovial has printed graduations and feature a white marking area for sample identification. Caps and tubes have both the same amount of expansion, which further enhances their leakproof qualities at changing temperatures.

The caps are manufactured in polyethylene and embody a silicone O-ring that guarantees them to be leakproof.

Colour code inserts can be fitted into the caps for identification.


All cryovials are certified RNase, DNase, and pyrogen free.

They are sterilized by radiation, packaged in safety-lock bags of 100, and are supplied capped.

They can withstand centrifugation up to 14,000xg.

Specifications:

Volume (ml)	2
Skirt	No
Dimensions (mm)	12.5 x 48
Material	Polypropylene
Quantity	1000

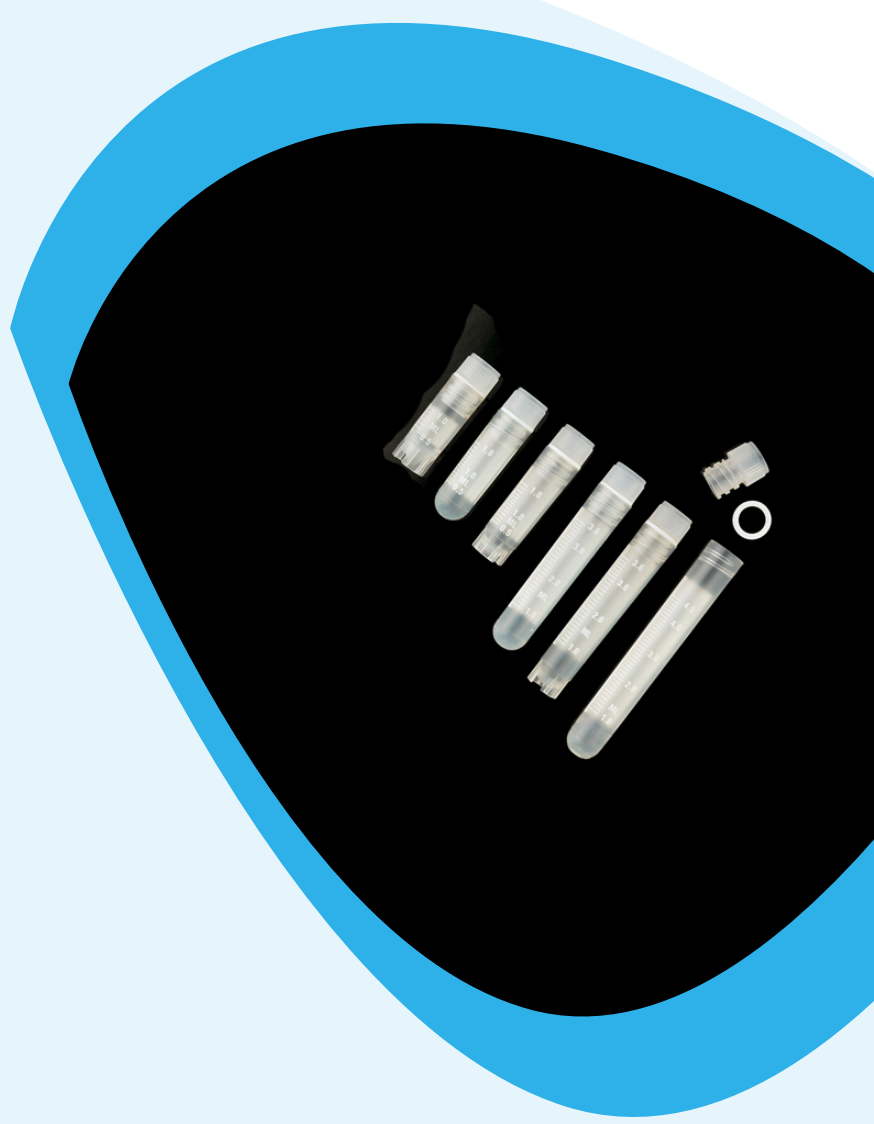

Unit 3, Tower Business Park
Warpsgrove Lane
Chalgrove
Oxfordshire
OX44 7XZ

Tel: (+44) 01865 400321

Email: enquiries@medlinescientific.com

Website: www.medlinescientific.com