

## Deltalab 1.5ml Premium Microtubes with Safety Lock

**Product Code:** 4092.5N

**Brand:** Deltalab

Made of ultra clear polypropylene, Deltalab's 1.5ml Premium Microtubes with Safety Lock feature a special "snap-seal" secure lock.


A special locking clasp with a distinctive audible snap-seal design ensures the tubes will not open during centrifugation. Each microtube has a frosted writing surface on the side and on the attached closure for better sample identification.

The caps are easy to handle and can be opened and closed with one hand. The lids are easily pierceable. The tubes have moulded-in graduations in 0.5ml increments.

They are autoclavable and can be used from -90°C to 121°C. They can withstand centrifugation up to 20,000xg and are free of DNA, DNase, RNase, and pyrogens.

### Specifications

<b>Dimensions (Ø x h; mm)</b>	11 x 39
<b>Quantity</b>	500

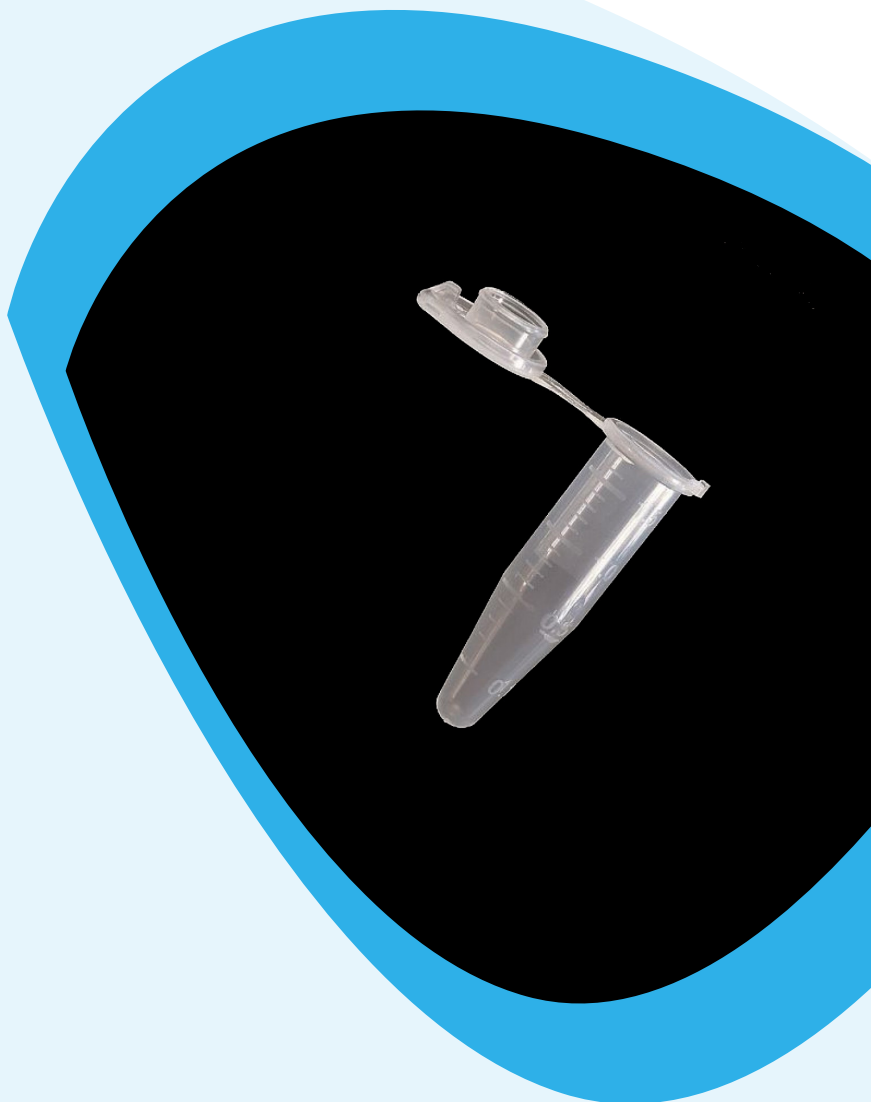

Unit 3, Tower Business Park  
Warpsgrove Lane  
Chalgrove  
Oxfordshire  
OX44 7XZ

**Tel:** (+44) 01865 400321

**Email:** [enquiries@medlinescientific.com](mailto:enquiries@medlinescientific.com)

**Website:** [www.medlinescientific.com](http://www.medlinescientific.com)