

Deltalab White Tissue Embedding Cassettes for Biopsy (Histosette I)

Product Code: 554111


Brand: Deltalab

Made of polyacetal, Deltalab's White Tissue Embedding Cassettes for Biopsy (Histosette I) have air vents allowing more efficient filling with paraffin. They have a writing area at a 45° angle and an incorporated one-piece integral lid.

Each cassette has 224 draining holes (1 x 1mm). They also have a safety closure. To lock the cassette, just snap apart the lid, turn it round, and lock it into the base of the cassette. They may be opened and closed as many times as required.

Specifications

External Dimensions (mm)	40 x 28 x 6 (with closed lid)
Colour	White
Quantity	1500

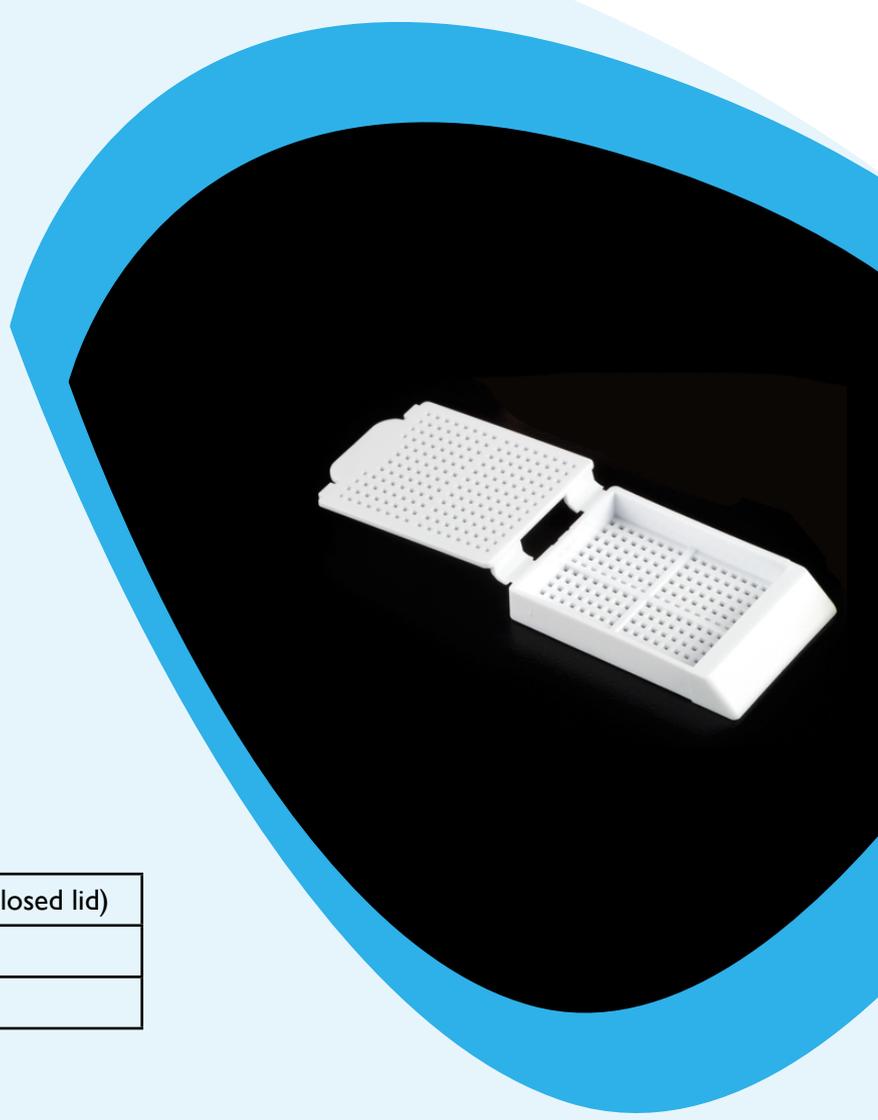

Unit 3, Tower Business Park
Warpsgrove Lane
Chalgrove
Oxfordshire
OX44 7XZ

Tel: (+44) 01865 400321

Email: enquiries@medlinescientific.com

Website: www.medlinescientific.com