

## Deltalab Assorted CryoInstant

**Product Code:** 409113/6

**Brand:** Deltalab

Deltalab's Assorted CryoInstant is an excellent preservation system for storing microbiological culture (e.g. fungus or sporulation phase), consisting of 2ml cryovials containing 25 porous beads and cryopreservative-added broth.

With these vials, we can:


- Get a perfect freezing storage system
- Obtain up to 25 replicates from the original microorganism (to use them for many years)
- Ease the recovery of the stored strain streaking an incubated bead directly onto solid media
- Prevent the changes and repetitive thawing/freezing of the rest of the beads
- Avoid the growth of ice when recovering
- Improve recovery rates of frozen cultures by adding cryopreservative to broth
- Minimise the risk of cross contamination
- Save space in the freezers when storing.

The external thread cryovials are made of polypropylene. Their smooth inner surface prevents from contamination. Their caps embody a silicone washer to ensure a positive leakproof seal. Tubes and caps withstand up to  $-190^{\circ}\text{C}$ .

Cryovials are presented in carton boxes (resistant to  $-100^{\circ}\text{C}$ ; dimensions: 150 x 150 x 55mm) with 100 units. Each box is supplied plastic-wrapped and labelled with code reference, batch, expiration date, and specifying the colour of the caps and beads.

### Specifications

<b>Cap &amp; Bead Colour</b>	Assorted - 5 colours x 20 cryovials
<b>Quantity</b>	100

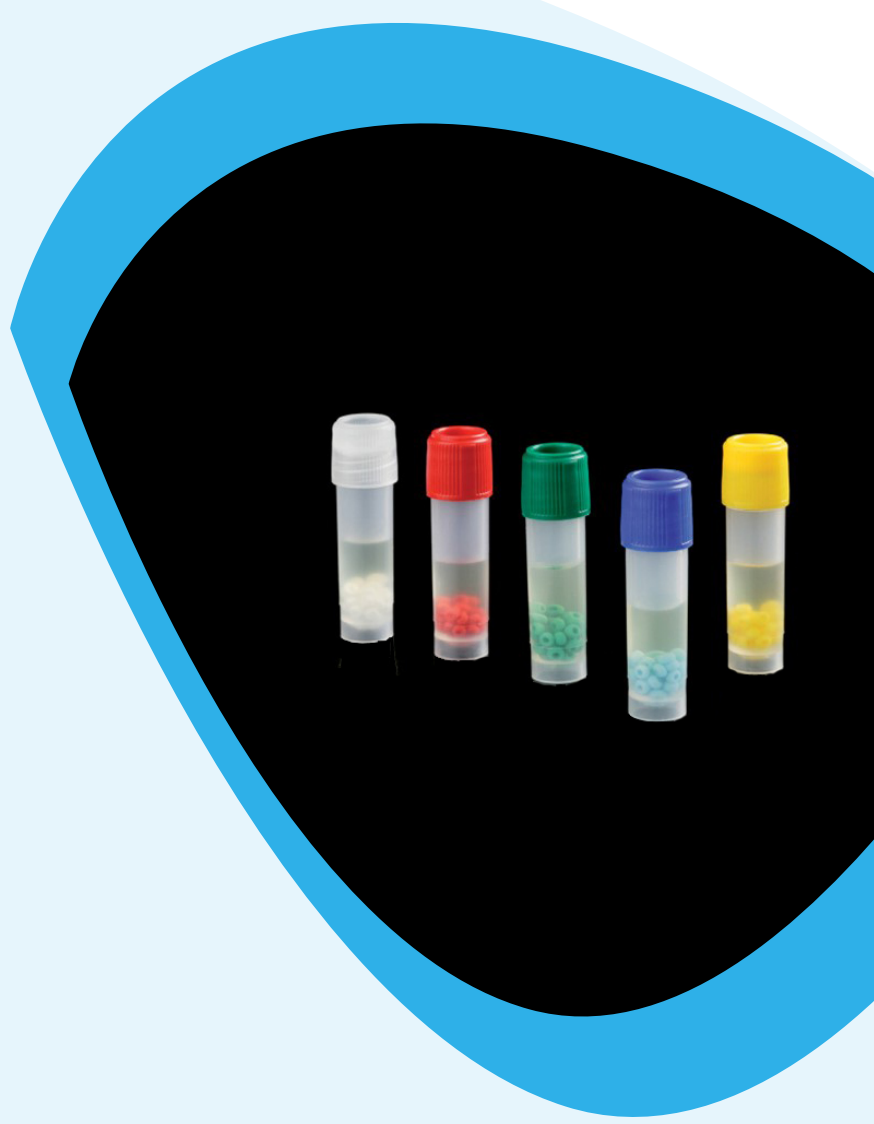

Unit 3, Tower Business Park  
Warpsgrove Lane  
Chalgrove  
Oxfordshire  
OX44 7XZ

**Tel:** (+44) 01865 400321

**Email:** [enquiries@medlinescientific.com](mailto:enquiries@medlinescientific.com)

**Website:** [www.medlinescientific.com](http://www.medlinescientific.com)