


## Jeio Tech 8 to 100rpm Waving Shaker (CWS-350)

**Product Code:** AAH381115K

**Brand:** Jeio Tech

The Jeio Tech 8 to 100rpm Waving Shaker (CWS-350) comes with enhanced conveniences for tilt adjustment, removable tray, and digital control.

This laboratory shaker is ideal for ELISA, enzyme immunoassays, protein synthesis, hybridisation staining or destaining gels, pharmaceutical profiling, and other laboratory applications.

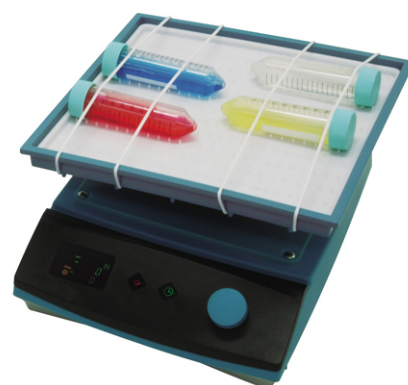

## Features

### Performance

- Microprocessor PID control
- Digital speed control providing a gentle to vigorous motion (8 to 100rpm)
- Easy to adjust tilt angle (0 to 13°)
- Smooth, quiet, uniform, and yet powerful shaking motion

### Convenience

- Vivid and colourful fluorescent display (VFD) with touch-sensitive keypad
- One-touch timer for time-dependent shaking operation (1 min. to 99 hr. 59 min.)
- Compact and lightweight design allow for use in cold rooms or incubators
- Three types of removable trays for easy transport of well-plates, test tubes, etc. (standard dimpled mat, and optional rubber mat and universal harness)
- Variety of accessories are available to meet application demands  
Accommodates to be used in a variety of test tube size, beakers, vials, micro well-plates, culture dishes, and plastic/glass trays

### Safety

- Self-diagnostic function identifying errors
- Its polypropylene main body is highly chemical-resistant and easy-to-clean
- Over-current protection

## Specifications

<b>Control System</b>	PID feedback control
<b>Display</b>	VFD
<b>Motion</b>	Waving
<b>Speed Range (rpm)</b>	8 to 100
<b>Tilt Angle</b>	1 min. to 99 hr. 59 min.
<b>Max. Load (kg)</b>	Up to 5.2
<b>Tray Dimensions (W x D)</b>	276 x 276mm
<b>Dimensions (W x D x H)</b>	249 x 325 x 138mm
<b>Weight (kg)</b>	3.7

Unit 3, Tower Business Park  
Warpsgrove Lane  
Chalgrove  
Oxfordshire  
OX44 7XZ

**Tel:** (+44) 01865 400321

**Email:** [enquiries@medlinescientific.com](mailto:enquiries@medlinescientific.com)

**Website:** [www.medlinescientific.com](http://www.medlinescientific.com)