


Deltalab Absorbent Paper Reel

Product Code: I20003

Brand: Deltalab

Deltalab's Absorbent Paper Reel provides total protection thanks to its two layers. The top side consists of a layer of cellulose with great liquid absorption power; the reverse side consists of a layer of polyethylene ensuring that it is completely waterproof. This absorbent paper can be used for the following applications: laboratory trays and tables protection, pathological anatomy laboratories, seed germination (between-paper method), and spillage recovering using its leaktight surface. It is particularly suitable for working with materials that are valuable (precious metals) or hazardous (toxic or corrosive materials).

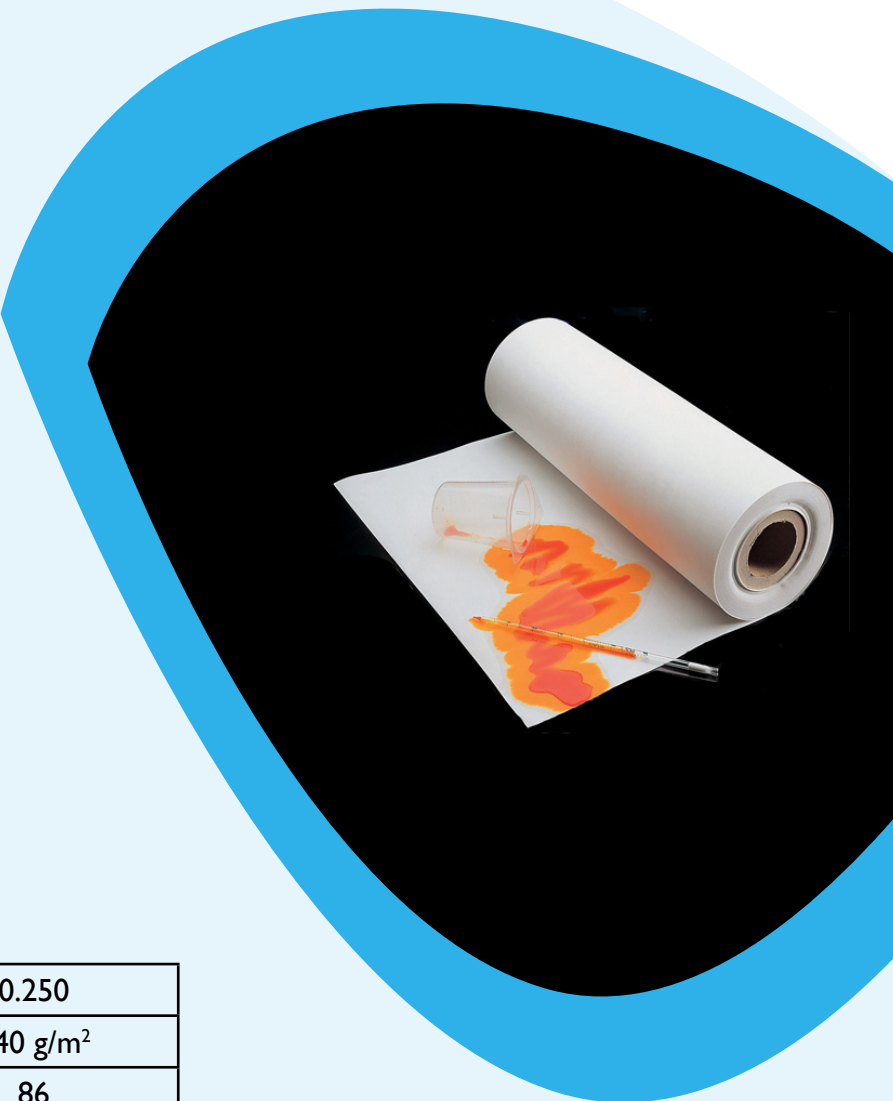

Specifications:

Paper Thickness (mm)	0.250
Absorption	240 g/m ²
Klemm Absorption (mm)	86
Description	1 reel
Material	Paper
Sterility	Non-Sterile
Quantity	1

Unit 3, Tower Business Park
Warpsgrove Lane
Chalgrove
Oxfordshire
OX44 7XZ

Tel: (+44) 01865 400321

Email: enquiries@medlinescientific.com

Website: www.medlinescientific.com