

Jeio Tech 300 Litre Plant Growth Chamber (GC-300TLH)

Product Code: AAHA1022K

Brand: Jeio Tech

Gain guaranteed photosynthesis conditions with the Jeio Tech 300 Litre Plant Growth Chamber (GC-300TLH).

With a capacity of 300 litres, this plant growth chamber delivers an innovative performance of low and high temperature along with humidity control. This is attained by altering the humidification tank exterior to the chamber, giving a joined advantage for carrying out frequent tests.

Features:

Performance:

General Control System

- 5°C to 50°C (lamp off)/10°C to 50°C (lamp on).
- Max. to 20,000 Lux.
- 40 to 80% RH (at 20 to 35°C).
- Max. 5000ppm CO₂ on/off system (optional).
- Microprocessor PID control/temperature calibration/automatic tuning.
- 10 step programmable temperature, humidity, and illumination profiles and repeatable steps of up to 999 cycles.
- Advanced performance of low and high temperature and humidity control, achieved by adapting the humidification tank exterior to the chamber; added benefit for long term test.

Illuminance Control System

- Uniformed luminescence distribution.
- Broad distribution of side lamps for hastening the growth of plants.
- Unique construction for minimization of heat increase from surrounding lamps. Tempered glass door blocks heated air from lamps.
- Stable and long lasting lighting through introduction of high frequency electronic ballast lamps.

Convenience:

- Ergonomically designed door construction. Inner glass door with silicone and external door with magnetic sealing for dual airtight packing. Well designed providing a smooth open for minimization of damage to plants.
- Wide inner tempered glass door for clear observation of plant growth without affecting inner chamber's environment.
- Tall plants can be grown inside of the chamber with adjustment of shelf level.
- Water supply during operation is also available for long term test purposes.
- Adjustable water tank position.
- Detachable condenser air-filter for easy maintenance of refrigerating efficiency. Maintenance of air filter no longer cumbersome with our detachable condenser air filter.
- Casters for easy mobility during installation or relocation.
- Eco-friendly CFC-free refrigerant use.


Safety:

- Automatic shut off after overheat alarm.
- Low and empty water level alarm.
- Power supply leakage breaker.
- Over current protection.
- Open door alarm.



Lamps Setting:

- Control illumination values with the below lamp setting arrangements.

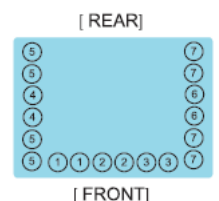

GC-300 lamp arrangement

Specifications:

Chamber Volume (litres)		300
Controller		Microprocessor PID controller
Temperature	Range (°C)	5 to 50 (lamp off)
		10 to 50 (lamp on)
		20 to 50 (with humidity)
	Fluctuation at 25°C	± 0.4 (without humidity)
	Variation at 25°C	± 1.0 (without humidity)
Illumination	Range (Lux)	0 to 20,000
	Control (steps)	10
	Lamp	FI lamp (32W x 18ea)
Humidity (RH%)	Range	40 to 80 (at 20 to 35°C) 70 to 90 (at 36 to 50°C)
	Fluctuation at 60% RH	± 3
Carbon Dioxide (optional)	Range (ppm)	300 to 5000
	Sensor	NDIR Carbon Dioxide
Temp. Heater Power (W)		1/2
Humid. Heater Power (W)		1800
Dimensions (w x d x h)	Internal	510 x 540 x 1100mm
	External	700 x 805 x 1900mm
Net Weight (kg)		250
Electrical Requirements (230V)		50Hz, 15A

Unit 3, Tower Business Park
Warpsgrove Lane
Chalgrove
Oxfordshire
OX44 7XZ

Tel: (+44) 01865 400321

Email: enquiries@medlinescientific.com

Website: www.medlinescientific.com