


## Medline Mini Vortex Mixer (FE101001)

**Product Code:** FE101001

**Brand:** Medline Scientific

The small but effective Medline Mini Vortex Mixer (FE101001) is a robust unit, specially designed for mixing tubes. To start mixing, the tubes just need to be pushed gently into the mixer's soft attachment. Mixing stops as soon as the tube is removed from the attachment.

With its unique foot pad design, this vortex mixer can even mix full 50ml containers without "walking". Its compact design makes it easy to store and use the mixer in small spaces (e.g. fume hoods). This would make it ideal for personal use and for use in teaching laboratories. It can be used in cold rooms and incubators (non CO<sub>2</sub>).

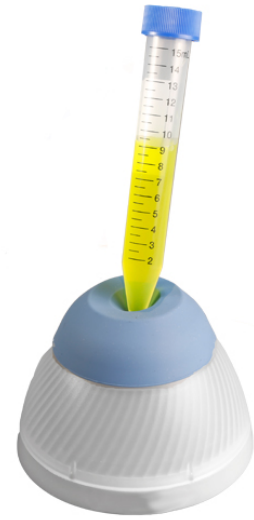

### Features

- Touch function
- Fixed speed of 3000rpm
- Orbital diameter of 5.5mm
- Excellent mixing capability can make full 50ml containers mix within 5 seconds
- Patented foot design can hold the vortex mixer firmly in place when mixing
- Low operation voltage (12V) for full protection
- Power supply is CE, cULus, CCC, SAA, and FCC certified

## Specifications

<b>Shaking Motion</b>	Orbital
<b>Orbital Diameter (mm)</b>	5.5
<b>Speed (rpm)</b>	3000 (fixed)
<b>Control</b>	Tube Pressure
<b>Dimensions (Ø x h)</b>	100 x 78mm
<b>Weight (kg)</b>	0.55

Unit 3, Tower Business Park  
Warpsgrove Lane  
Chalgrove  
Oxfordshire  
OX44 7XZ

**Tel:** (+44) 01865 400321

**Email:** [enquiries@medlinescientific.com](mailto:enquiries@medlinescientific.com)

**Website:** [www.medlinescientific.com](http://www.medlinescientific.com)